



Provide Warmth  
Anywhere, Everywhere

Embrace Nest is an easy-to-use, portable infant warmer that does not require continuous power and is used as supplement to KMC.

NICU



Transport



Maternity Ward





# Component Features

## Embrace Nest Unique Benefits




**WARMPAK™**

Maintains temperature of ~37°C for at least 4 hours\*

\*duration of at least 4 hours at ambient temperature greater than 25°C

380 mm x 220 mm • 1.3kg

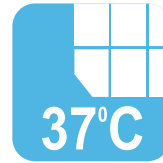

Special Phase Change Material in WarmPak™ maintains a temperature of ~37° C for at least 4 hours\*



Portable and simple to use



**BABYWRAP™**

Reusable and easy to sanitize using disinfectants or soap wash

520 mm x 250 mm • 0.3kg



Only 15 Watts of power consumption in ready mode



**ACCUTEMP HEATER**

Precision heating mechanism prevents overheating of WarmPak™

440mm x 290 mm x 70 mm • 2.5kg • 230 V, 50 Hz




Does not require a constant supply of electricity



Easy to sanitize and reuse

## Simple To Use

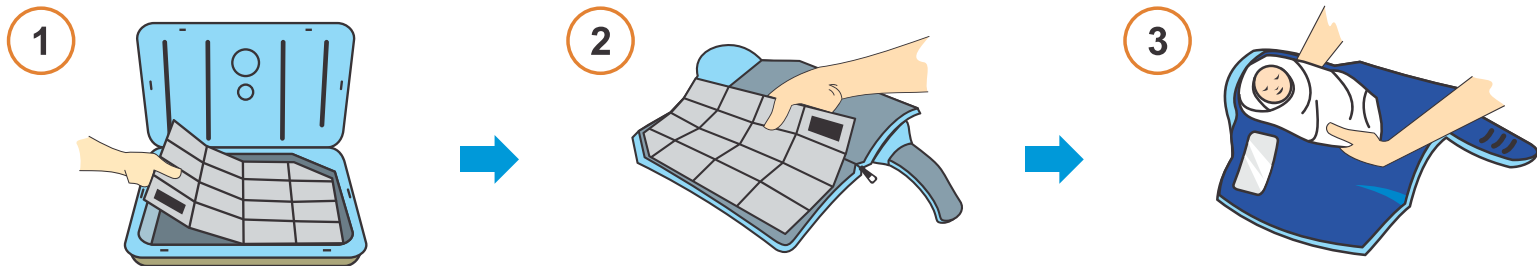

## Product Order Information

Model no. (FD No.)	Product Name
FD0904	Embrace Nest Starter Pack
FD0903	Warmpak™ Refill Pack
FD0902	Accutemp Heater Refill Pack
FD0901	BabyWrap™ Refill Pack

## Phoenix Medical Systems Private Limited

DP42, SIDCO Industrial Estate, Thirumudivakkam, Tamil Nadu – 600132, India

Toll Free India : 1800 425 3503, Phone: +91 44 2478 1705 / 2227,

✉ India : sales@pmsind.com, International: exports@pmsind.com

🌐 www.phoenixmedicalsystems.com